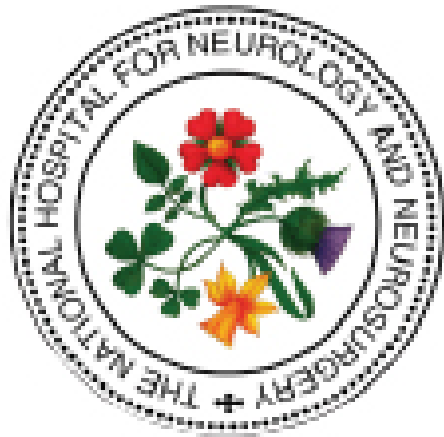




The National Hospital for Neurology and Neurosurgery Neuroscience Critical Care Unit (NCCU) LIQUOGUARD AUTOMATED INTRACRANIAL PRESSURE CONTROLLED CSF DRAINAGE

1. The drain is set up by the neurosurgeon

Ensure **red bung** replaces vented bung

2. Keep plugged in and do not switch LigoGuard machine off when in use at bedside

- If LigoGuard loses power or is switched off, it will return to factory settings and for some patients they may drain too much CSF – **machine will function on battery power for transferring patient**

3. Zero reference point

- As with all external ventricular drainage systems **zero reference point** is **tragus**
- With LigoGuard system place **ECG dot** at temple and attach sensor to set zero reference point

4. Settings

- Drainage mL/h and pressure mmHg
- Will only drain when actual pressure (**white curve**) is higher than target pressure (**blue line**) – therefore it is possible CSF drainage will be less than set flow rate

6. Monitor drainage volume (record hourly)

- Select **Info** in bottom tabs and then select **volumes** from right hand menu
- The **pSet** (Pressure set) and **vSet** (Volume set) should be checked and documented every hour along with the **volume drained**

5. Alarms - If alarm goes off

- Silence / Troubleshoot / Seek help
- Pulsation alarms can be disabled when used as a **lumbar drain** as the waveform is flatter and the LigoGuard machine may keep alarming

High/Low	May trigger if patient sits up too rapidly – raise head of bed slowly
Pulse	Check if blocked or kinked
Flow	Too much / too little drainage

7. Taking patient to CT scan

- Place in pause mode or keep running dependent on individual patient

8. Raising or lowering height of bed / sitting patient out of bed

- No need to pause – drainage is not dependent on relation of machine to zero reference point

9. Special considerations

- If pump alarms to say drainage set has expired (*manufacturer's advice change every 7 days*) but neurosurgeon happy to keep same set attached (usually up to 10 days), try changing LigoGuard pump but keep same tubing attached
- Empty drainage bag when full - no need to change entire set
- For MRI, drainage set is MRI compatible, LigoGuard machine and ECG electrode (used as sensor) are not** - Drainage set needs to be removed by medical staff from LigoGuard pump – replace after MRI without need to disconnect drain - *not required for CT / xray*

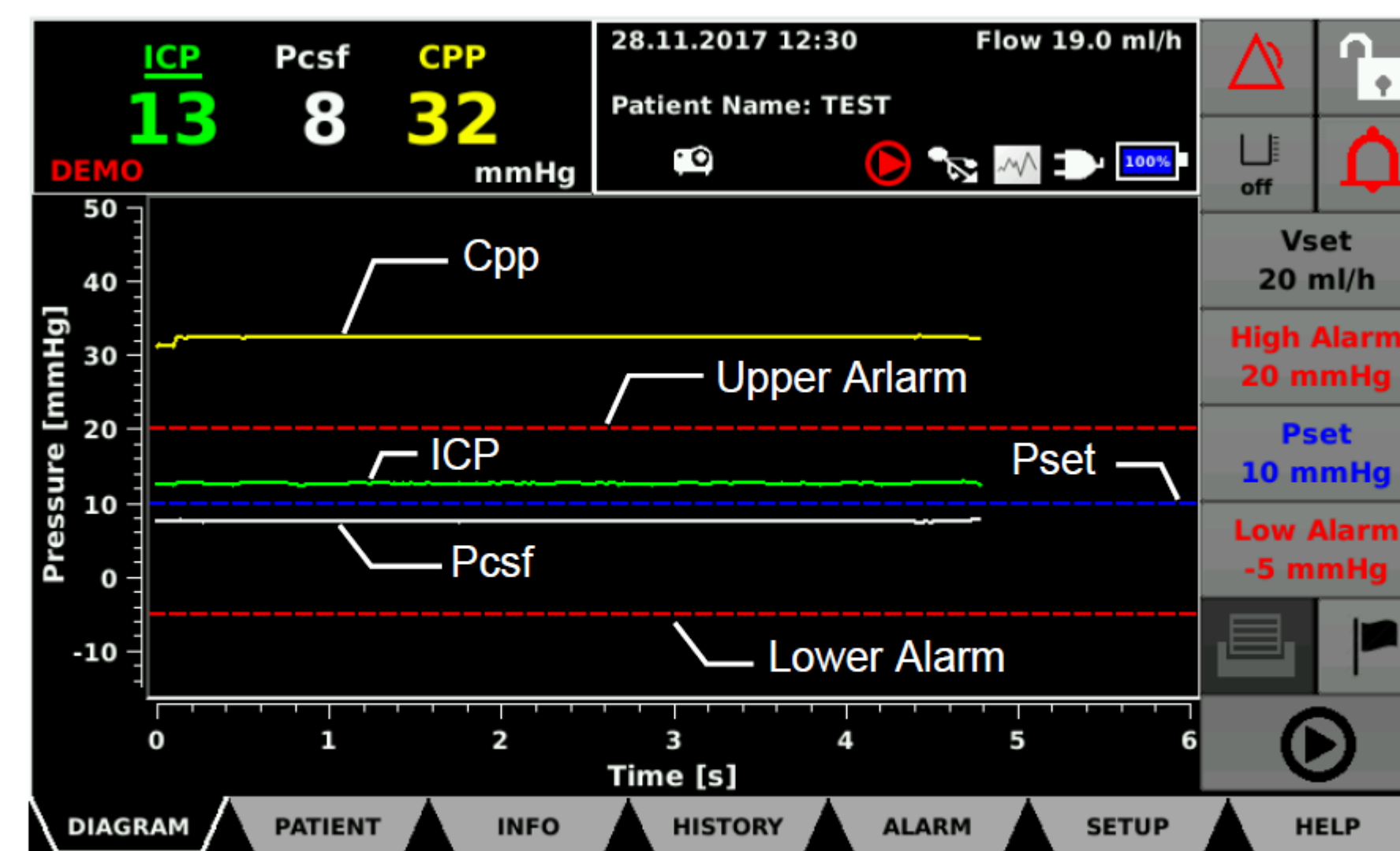
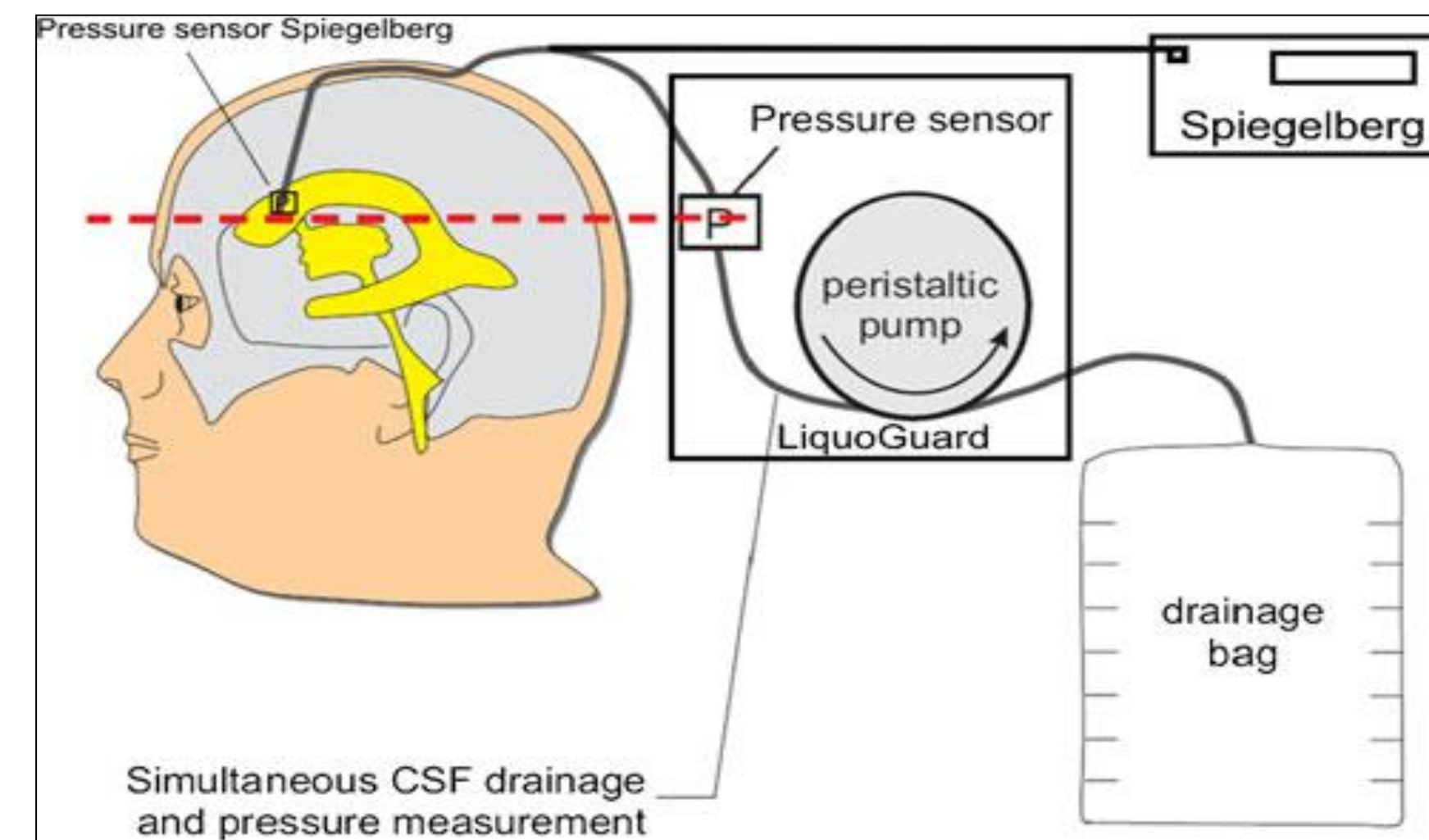


Figure 11: Limit values and curves shown in the diagram.

- Pressure (y-axis) in mmHg or cmH₂O
- White curve: CSF pressure measured by pressure sensor in tubing set
-
- Green curve: ICP pressure value measured by an external sensor (parenchymal or catheter tip = ventricular)
- Yellow curve: cerebral perfusion pressure (CPP)
- Blue dotted line: **Pset** When CSF pressure (ICP or Pcsf) higher than Pset, and application not paused, pump turns (drainage) to reduce CSF pressure to Pset
- Top red dotted line: **Upper Alarm** (upper pressure alarm limit)
- Lower red dotted line: **Lower Alarm** (lower pressure alarm limit)

CSF SAMPLING

LigoGuard set does not have a drip chamber therefore CSF samples must be taken from the 3-way tap – **follow CSF sampling protocol**

ADMINISTERING IT ANTIBIOTICS

LigoGuard machine must be set to **PAUSE MODE** when administering IT antibiotics – **follow IT protocol**

